


Fit-N-Swipe

Cleaning Pads with Purpose



- High cleaning effect and highly absorbent
- Practical, single-handed use
- Can be positioned individually in the direct working area
- Fit-N-Swipe Clean ideal for instruments
- Fit-N-Swipe Dry ideal for mirrors

Can be placed anywhere!

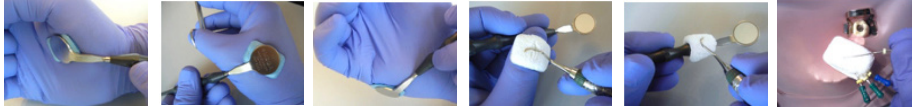

Ball of the thumb The hand Heel of the hand The finger The instrument The rubber dam


Fit-N-Swipe Clean

Adhesive cleaning pad for clean instruments



- Unique microfiber with high retention capacity
- Safe and handy pad, use as pincushion for endodontic needles

605251	Fit-N-Swipe, Clean (white)	50 pads
--------	----------------------------	---------

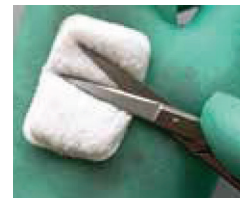

Cleaning scissors




Cleaning scaler



Use as a pincushion for endodontic needles


Fit-N-Swipe Dry

Adhesive drying pad for crystal clear mouth mirrors

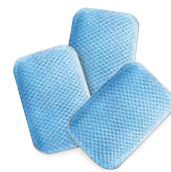


- Unique multi-layered fabric
- Removes even greasy or oily contaminants without streaking

605252	Fit-N-Swipe, Dry (blue)	50 pads
--------	-------------------------	---------

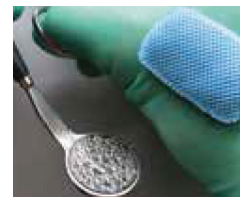

Wet and fogged mouth mirror

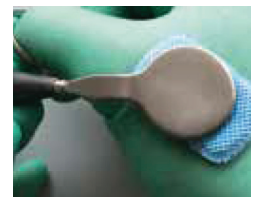

Simply wipe off on the pad



Crystal clear mouth mirror

Testimonials

Bettina Drewes, DH: "I think Fit-N-Swipe is a very good product as there usually is no assistant available in prophylactic treatment and Fit-N-Swipe makes the work so much easier."

Melanie Graf, ZMP: "Absolutely great! It's simple and easy. 'Peel off and stick on the glove - ready!' After the treatment simply dispose of the glove."

Check out the video



* Can be sterilized in autoclave before use



Fit-N-Swipe

Instructions for use

Fit-N-Swipe Clean

Adhesive pads to help clean instruments during use

- Convenient single-handed use
- Unique microfiber fabric with high retention capacity
- High cleaning effect and highly absorbent
- Versatile placement options
- Pads are single-use and non-sterile. They can be sterilized prior to the single-use in an autoclave according to the instructions for use for textile fabrics.

Positioning the Pads

- Bend the pad for an easier removal of silicone paper
- Press the pad firmly. When changing position a new pad must be used
- All surfaces (glove, rubber dam) must be absolutely dry
- Any residue of disinfectant or soap must be removed
- Do not stick the pad directly to the skin and do not reuse it

Adhesion

Adhesion tested on vinyl, nitrile, latex. Special glove coatings can reduce the adhesiveness.



Bend the pad for an easier removal of silicone paper

Fit-N-Swipe Dry

Adhesive drying pad for crystal clear mouth mirrors

- Unique multi-layer fabric
- High cleaning effect and highly absorbent
- Dries mirror streak-free
- Versatile placement options
- Pads are single-use and non-sterile. They can be sterilized prior to the single-use in an autoclave according to the instructions for use for textile fabrics.

Positioning the Pads

- Press the pad firmly. When changing position a new pad must be used
- All surfaces (glove, rubber dam) must be absolutely dry
- Any residue of disinfectant or soap must be removed
- Do not stick the pad directly to the skin and do not reuse it

Adhesion

Adhesion tested on vinyl, nitrile, latex. Special glove coatings can reduce the adhesiveness.



Bend the pad for an easier removal of silicone paper

